

Invia® Liberty™ NPWT System

THERAPY CONFIDENCE MADE SIMPLE

Therapy Confidence

- Intelligent Pressure Control from the pump to the wound site ¹
- Dynamic Exudate Removal helps prevent blockages from the wound site to the pump
- Secure canister and tubing connection² with an audible click

Simplified Management

- Intuitive pump operation and troubleshooting instructions
- Compact design for patient comfort and mobility
- Clinical flexibility with comprehensive dressings portfolio



Medela provides efficient solutions to support clinicians with Negative Pressure Wound Therapy management in all patient care settings.

Invia Liberty delivers therapy confidence by intelligently controlling set pressure at the wound site 1. The pump dynamically adapts the air flow cycle to both volume and viscosity of fluid in the tubing to optimize exudate removal.

Click-in canister attachment and one-touch canister release button allow for easy canister exchange. Secure tubing attachment is guaranteed with the new Invia Quick-connector feature².






Intuitive and quiet pump operation promotes full integration in the patient's daily life, while patient comfort and mobility are facilitated by Invia Liberty's compact design and convenient carrying case.

Key features at a glance

- Intelligent Pressure Control from the pump to the wound site
- Dynamic Exudate Removal actively helps prevent blockages with air flow cycles adapting up to four times a minute
- Pressure options from -40 to -200 mmHg including standard -125 mmHg
- On-screen air leak indicator
- Custom air leak modes for normal and large wounds
- Easy and secure tubing attachment with the Quick-connector
- One-touch canister release button
- Choice of 300 ml and 800 ml disposable canisters with solidifier



Technical data

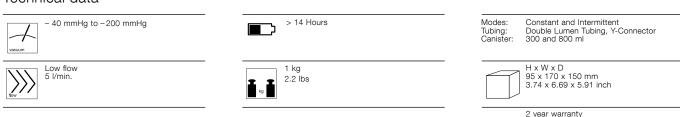

Find out more about the new Invia Liberty system at www.medelahealthcare.com

1	Performance data on file at Medela AG
2	Usability Validation data on file at Medela

Medical Vacuum Technology for Healthcare Professionals

Please contact us or your local Medela representative for details.

ocal contact:			

Medela, Invia and Invia FitPad are trademarks of Medela Holding AG, Switzerland.



United Kingdom Medela UK Ltd. Huntsman Drive Northbank Industrial Park Irlam, Manchester M44 5EG United Kingdom Phone +44 161 776 0400 +44 161 776 0444 info@medela.co.uk www.medelahealthcare.co.uk International Sales Medela AG Lättichstrasse 4b 6341 Baar Switzerland Phone +41 41 562 51 51 +41 41 562 51 00 customercare@medela.ch www.medelahealthcare.com